

M-34 Digital PP Rev.C

microson



- ① Ear Hook
- ② Microphone
- ③ Switch (M/T)
- ④ Volume Control
- ⑤ Battery compartment (Type 13)
- ⑥ Trimmers cover (NH, PC)

DESCRIPTION

Behind The Ear air conduction hearing aid (BTE), 13 Battery, M-T Switch (Microphone-Telecoil) 2 trimmers (NH,PC) and Volume Control. Battery Door with ON/OFF Functionality.

INTENDED USE

The MICROSON M-34 Digital PP (Rev.C) hearing instrument is indicated to compensate transmissive severe hearing losses. It is not suitable for children.

See Fitting Range⁽¹⁾


Features

- ✓ Matrix 136/60 @ 2cc
- ✓ 100% Digital Technology
- ✓ 2 Channels
- ✓ 2 Trimmer Controls (NH, PC)
- ✓ Switch (M/T)
- ✓ Volume Control
- ✓ Telecoil
- ✓ Feedback canceller
- ✓ Battery compartment with ON / OFF function
- ✓ Low Battery Indicator
- ✓ Small Size
- ✓ Fully compatible with mobile phones*
- ✓ Battery 13 Type - PR48 (IEC 60086)

Requirements


- ✓ No programming unit required

Weight & Dimensions

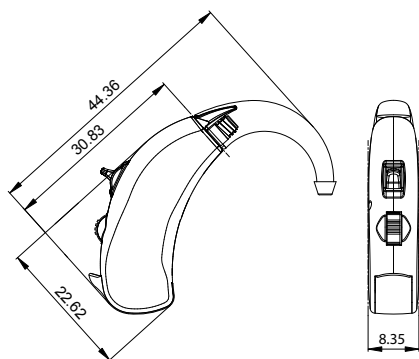

Weight excluding battery: 3.6 gr.
Weight including battery: 4.44 gr
Dimensions in millimeters (mm)

Fitting Range⁽¹⁾

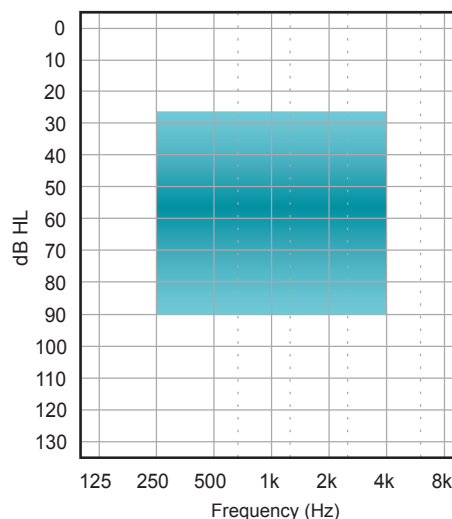

*In compliance with IEC 60118-13

DS-0012-020-EN
Rev.D 2017/01/09

M-34 Digital PP Rev.C
1 of 4

	Acoustic Data	IEC 60118-7:2005	IEC 60118-0:1993/A1:1994
OUTPUT	OSPL ¹ 90 Peak (dB _{SPL})	136	139
	OSPL90 Peak Frequency (Hz)	1200	1300
	HFA ² -OSPL90 / RTF ³ -OSPL90 (dB _{SPL})	129	135
GAIN	HFA-FOG ⁴ (dB)	52	60
	RTF-FOG (dB)	50	59
	FOG (dB)	60	64
	FOG Frequency (Hz)	1200	1300
	RTG ⁵ (dB)	51	59
NOISE	Equivalent Input Noise (dB _{SPL})	21	20
AGC ⁶	Attack Time (ms)	1	2
	Release Time(ms)	17	5
TELECOIL	HFA-SPLI ⁷ / RTF-SPLI (dB _{SPL})	111	119
	HFA-FOG-MASL ⁸ / RTF-FOG-MASL @ 1 mA/m (dB _{SPL})	83	89
	500 Hz @ 100 mA/m (THD %)	N/A	N/A
	800 Hz @ 100 mA/m (THD %)	2,8	6,6
	1600 Hz @ 100 mA/m (THD %)	2,4	3,7
DISTORTION	500 Hz @ 70 dB _{SPL} (THD %)	1,2	1,8
	800 Hz @ 70 dB _{SPL} (THD %)	0,7	1,4
	1600 Hz @ 65 / 70 dB _{SPL} (THD %)	0,2	0,2
CONSUMPTION	Current Drain (mA)	0,64	0,60
FREQUENCY LIMITS	f ₁ (Hz)	100	200*
	f ₂ (Hz)	4900	4900*

Power Source: 1.3V Battery Simulator

IEC 60318-5:2006 IEC 60318-4:2010

¹OSPL: Output Sound Pressure Level

²HFA: High Frequency Average

³RTF: Reference Test Frequency (1600 Hz)

⁴FOG: Full On Gain

⁵RTG: Reference Test Gain


⁶AGC: Automatic Gain Control

⁷SPLI: Sound Pressure Level Inductive


⁸MASL: Magneto Acoustical Sensitivity Level

* According to DIN 45605 standard


OSPL90 / OSPL50 @ FOG @ IEC 60118-7:2005

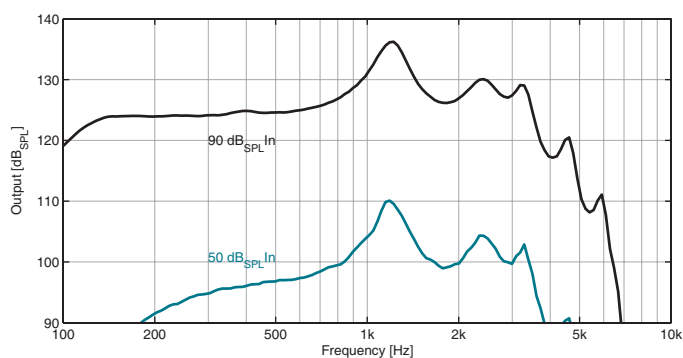


OSPL90 / OSPL50 @ FOG @ IEC 60118-0:1993/A1:1994

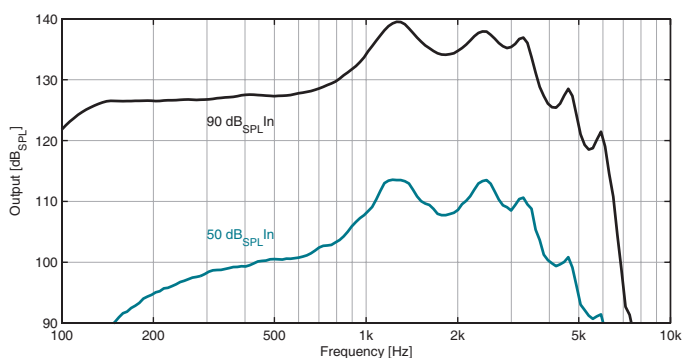

FREQUENCY RESPONSE @ RTS @ IEC 60118-7:2005

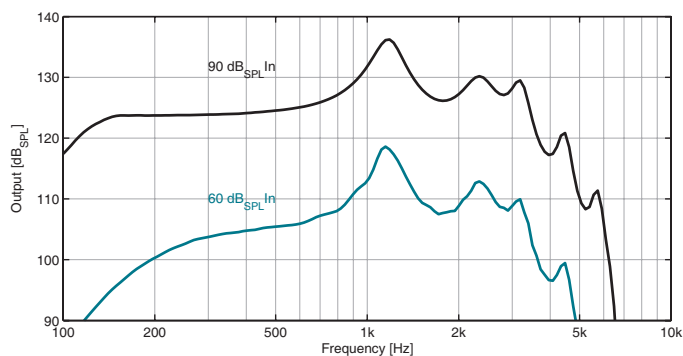


FREQUENCY RESPONSE @ RTS @ IEC 60118-0:1993/A1:1994

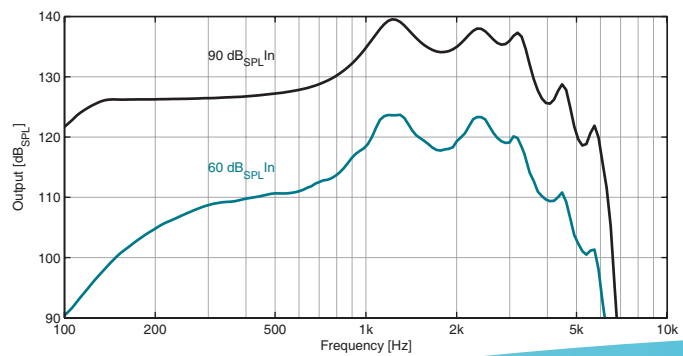
DS-0012-020-EN
Rev.D 2017/01/09


M-34 Digital PP

Rev.C
2 of 4


FREQUENCY RESPONSE @ VOL. CTRL. @ IEC 60118-0:1993/A1:1994

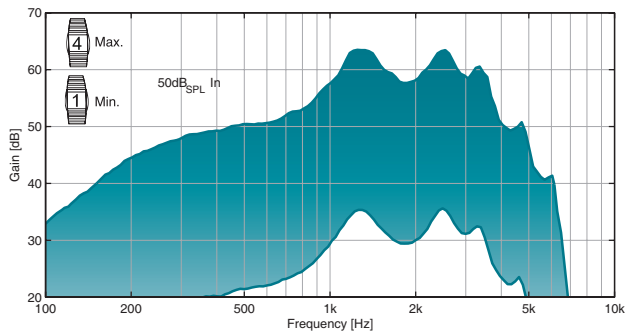


INDUCTION COIL SENSITIVITY @ RTS @ IEC 60118-7:2005

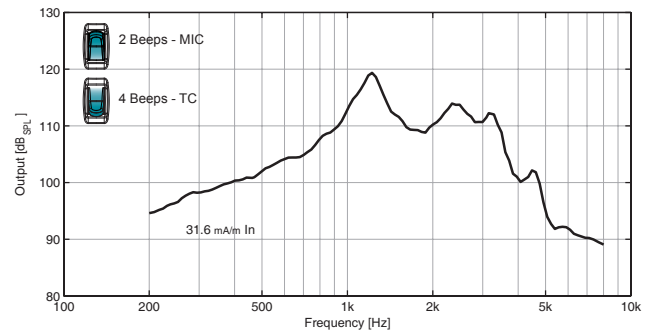


FREQUENCY RESPONSE @ NH CTRL. @ IEC 60118-0:1993/A1:1994

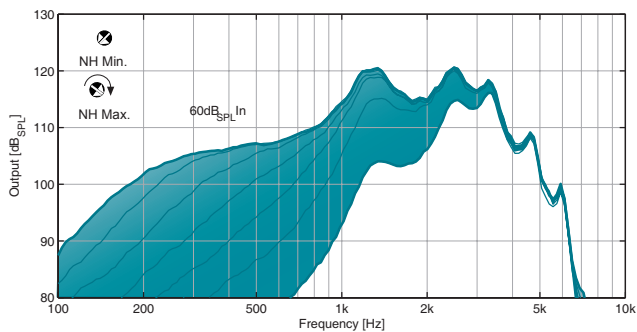


FREQUENCY RESPONSE @ PC CTRL. @ IEC 60118-0:1993/A1:1994

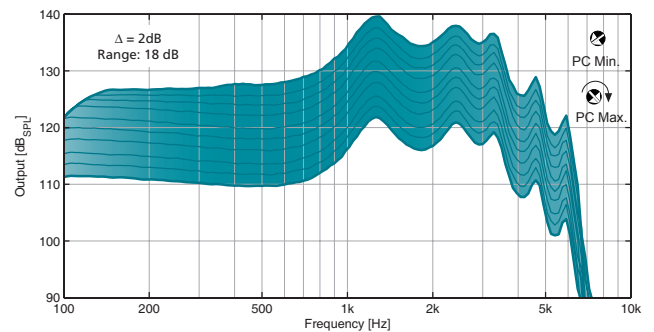

INPUT-OUTPUT @ RTS @ IEC 60118-7:2005

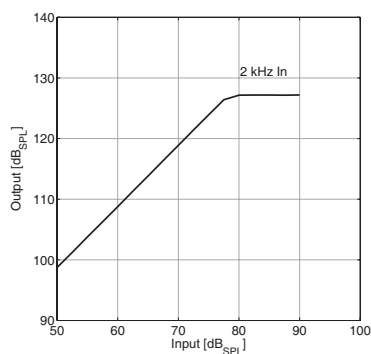

INPUT-OUTPUT @ RTS @ IEC 60118-0:1993/A1:1994

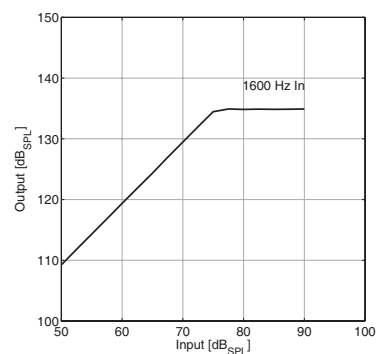

Spare Parts

91521 Earhook P13 M1 (5pcs)
 91212 M-34 Digital PP Blue&Red Label lid set
 74476 M-34 Digital PP Label Lid (10 pcs)
 79431 Windscreen BTE Beige (10 pcs)
 75140 Volume Control Cover (10 pcs)

User Manual

91360 User Manual M-34 Digital BTE (LP2) (ES/EN/IT/PT/FR)
 91588 User Manual M-34 Digital BTE (LP3) (EN/TR/RO/BG/EL/PL)

HEARING INSTRUMENT CLASSIFICATION IN COMPLIANCE WITH IEC 60601-1 STANDARD

Medical Device Classification

Protection against electric shock	MEDICAL DEVICE WITH INTERNAL ELECTRICAL POWER SOURCE
	B Type Applied Part
	 This symbol indicates that the products described in these user instructions adhere to the requirements for an application part of Type B of IEC 60601-2-66. The surface of the hearing aid is specified as an applied part of Type B.
Working Method	CONTINUED WORKING

Environmental conditions

	Temperature Min. (°C)	Temperature Max. (°C)	Relative Humidity Min (%)	Relative Humidity Max (%)
Recommended usage and storage	0	40	10	95

Power Supply Electrical Features

	M-34 Digital PP
Nominal Operating Voltage	1.4 V
Current Type	Direct current DC
Nominal Current Leakage	0,60 mA
Battery Nomenclature (IEC 60086)	PR48

PRODUCT	REFERENCE	MODEL	GTIN-13
M34 DIGITAL BTE	73230	M 34 DIGITAL PP	8435281304417

GMDN Code: 34671